

# 2°C~8°C Pharmacy Refrigerator

**5**-series High-quality Product

- Seven Sensors
- Temperature Uniformity
- Ensuring Safety of Storage



## YC-1015L



### Real Air Cooling

The multi-duct vortex refrigeration system and finned evaporator can prevent frosting completely and improve temperature uniformity to a large extent.

### Real Energy Saving

The refrigerator is equipped with “M + Energy Core” power management system and complete air cooling design, improving energy saving efficiency by more than 40%.

### Real Visibility

The electrical heating + LOW-E design with double considerations can achieve a better anti-condensation effect for the glass door.

### Real Safety

Seven temperature probes can ensure high precision of temperature control with nearly no fluctuation and thus can improve safety.





### Intelligent Temperature Control

- Equipped with high-precision microcomputer temperature control system and 1 high-brightness digital temperature display with display accuracy at 0.1°C.
- Equipped with a remote alarm access port and 485 data access port.




### Security System

- The well-developed audible & visual alarm system with multiple audible & visual alarm functions: high temperature alarm, low temperature alarm, sensor failure alarm, door ajar alarm, power failure alarm;
- The turn-on delay and stopping interval protection function can ensure reliability in running;
- The door is equipped with a lock, preventing it from opening accidentally.




### Refrigeration System

Equipped with a compressor supplied by an international famous brand, high-efficiency air-cooled condenser and finned evaporator, ensuring fast refrigeration.

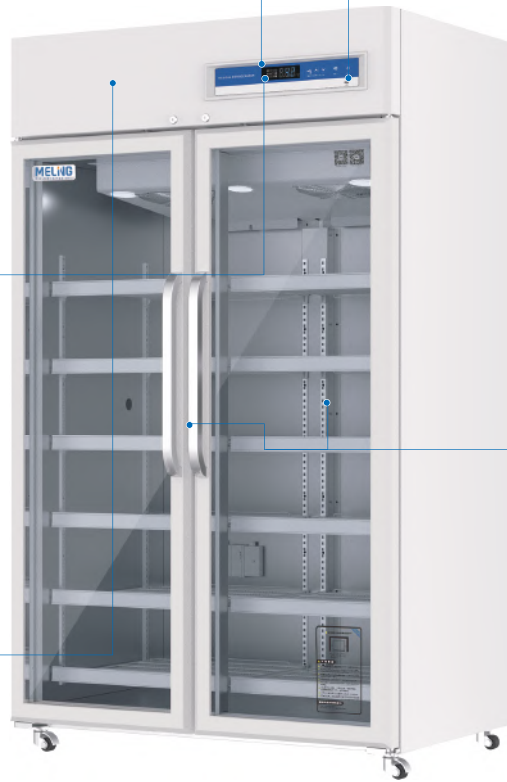

### Data Storage

- Equipped with a USB export interface, which can be used for exporting data in PDF;
- When a U-disk is connected, temperature data can be stored continuously and automatically (data can be stored for 100,000 sets of data);
- The printer is optional.

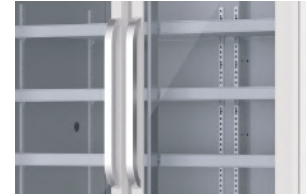


### Human-oriented

- Equipped with high-quality shelves made from PVC-coated steel wire with tag cards and classification marks, which are easy to clean;
- Upgrade with aluminium alloy door handle, ensuring elegance of appearance;
- Equipped with testing port, bringing convenience to users in testing temperature inside the cabinet.


## Scope of Application

Suitable for storage of biological products, vaccines, drugs, reagents, etc. Suitable for use in pharmacies, pharmaceutical factories, hospitals, centers for disease prevention & control, clinics, etc.

## Performance Data / Cooling Curve

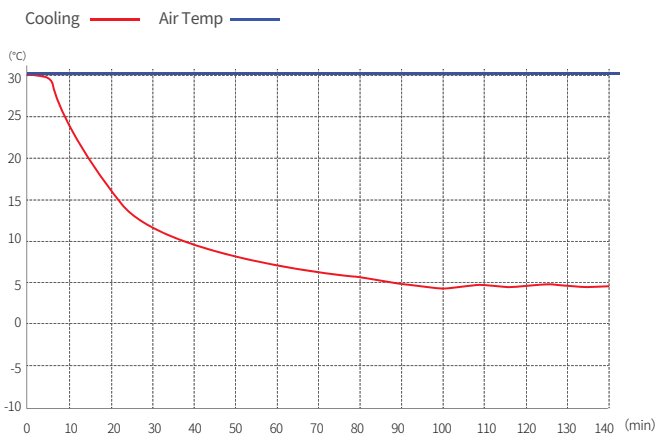

## External Dimensions

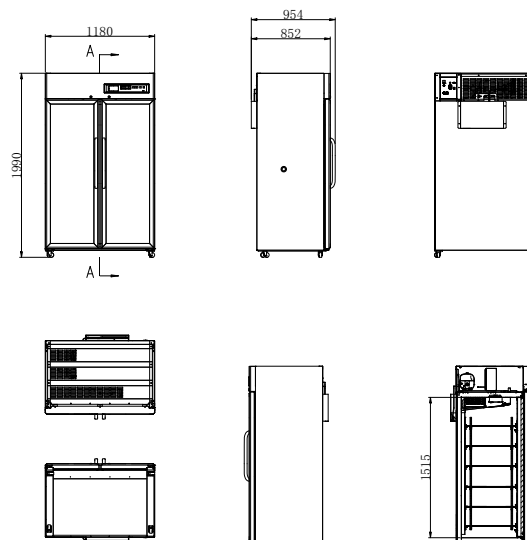

## Specification Chart

2~8°C Pharmacy Refrigerator	
Model	YC-1015L
Cabinet Type	Upright
Capacity(L)	1015
Internal Size(W*D*H)mm	1070*670*1515
External Size(W*D*H)mm	1180*954*1990
Package Size(W*D*H)mm	1313*988*2168
NW/GW(Kgs)	185/262
Performance	
Temperature Range	2~8°C
Ambient Temperature	16-32°C
Cooling Performance	5°C
Climate Class	N
Controller	Microprocessor
Display	Digital display
Refrigeration	
Compressor	1pc
Cooling Method	Air cooling
Defrost Mode	Automatic
Refrigerant	R290
Insulation Thickness(mm)	L/R/B:55,U:56,D:58
Construction	
External Material	PCM
Inner Material	Aluminum plate with spraying ( optional Stainless steel )
Shelves	12(coated steel wired shelf)
Door Lock with Key	Yes
Lighting	LED
Access Port	1pc. Ø 25 mm
Casters	4+ (2 with brakes)
Data Logging/Interval/Recording Time	USB/Record every 10 minutes / 2 years
Door with Heater	Yes
Alarm	
Temperature	High/Low temperature,High ambient temperature,Condenser overheating
Electrical	Power failure , Low battery
System	Sensor failure, Door ajar, Built-in USB datalogger failure, Communication failure
Electrical	
Power Supply(V/HZ)	220-230~/50
Rated Power(W)	520
Power Consumption(KWh/24h)	8.24
Rated Current(A)	3.21
Accessories	
Standard	RS485, Remote alarm contact, Backup battery
Options	Printer

\*The model, parameters and performance specified in this brochure may be changed without prior notice because of product upgrading.

\*There may be differences between the product images shown in this brochure and the actual products. When you are buying any product, please check the actual product.