

JIBI MED[®]

DEVOTE OURSELVES TO THE PROFESSIONAL MANUFACTURING OF THE STERILIZERS

| Company Profile

Jiangyin Binjiang Medical Equipment Co., Ltd.(JIBIMED) is a professional high-tech enterprise, incorporating research, manufacturing, distribution of sterilizing equipment and relevant medical products, located in Jiangyin city of Jiangsu Province, which is a remarkable place with picturesque scenery. Jiangyin is the birthplace of outstanding people, and ranked first on the economy among all the Chinese cities.

With the foundation concept of "technology reflects our wisdom, quality reflects our dignity". JIBIMED has gathered many top professional experts and famous consultants who involving themselves in the research and development of the sterile products after 10 years development. We have formed strong power to follow the tip technology at home and abroad, has become one of the largest and most professional manufacturers with continuous growth of export for the past years.

At present, JIBIMED has two producing zones, covering the area of 40 acres and the building area of 30000 m², we have more than ten full producing line for the vertical pressure steam sterilizer, portable pressure steam sterilizer, table top steam sterilizer, horizontal pressure steam sterilizer, low temperature plasma sterilizer, ultraviolet lamp trolley and intelligent drying oven, reaches 50000 sets for annual production, which are not only covered 70% of the Chinese hospitals, but also delivered to Europe, Africa, East-South Asia, Middle East over 50 countries, we have provided our goods and service to amount of tens of thousands of hospitals all over the world.

All the products from JIBIMED are strictly treated according to quality control system, such as CE, ISO9001, ISO13485. To supply the "zero-defect" units with comprehensive, considerate and sincere service is our final target goal.

The customers satisfaction and the quality expectation is our unswerving pursuit. Advanced technology, scientific management and excellent service is the guarantee of the continuously improving for the quality. Product safety, regulatory compliance and environmental protection is our social responsibility.

JIBIMED has win highly appreciation among domestic and oversea hospitals of various levels in the sterilizing field with innovated products, qualified service, and efficient management, has been awarded the "National Quality Trustworthy Enterprise", "Jiangsu Integrity and Outstanding Enterprises", ect .

For the honor and the commitment , we strive to do better !



Vertical Pressure Steam Sterilizer



Use

The vertical pressure steam sterilizers are assembled with heating system, micro-computer controlled system, over heat and over pressure protection system, which are reliable for the sterilizing effects, convenient of operation and energy conservation. They are suitable for sterilizing equipment, dressings, containers, liquids, culture medium, waste, etc. in hospitals, disease control centers, P2 biological laboratories and blood stations, etc. of scientific research institutes.

Characteristics

- ◆ Fully stainless steel structure.
- ◆ Hand wheel type of quick-open door structure.
- ◆ Door safety lock system.
- ◆ Automatically controlled by computer.
- ◆ LCD display of working status, touch type key.
- ◆ Over temperature & over pressure auto-protection.
- ◆ The internal steam/water circulation system is equipped with an exhaust cooling system, which avoids the generation of aerosols, and does not discharge steam to the outside. Therefore, the working environment is clean and dry.
- ◆ Self-inflating type seal.
- ◆ Safe protection of water lacking.
- ◆ Automatically shut off with beep reminding after sterilization.
- ◆ With two stainless steel sterilizing baskets.
- ◆ A high-efficiency filter is added to the exhaust to ensure the safety of the environment and operators.
- ◆ Optional drying system.

Fully automatic circulation
(specifically use for
clinical laboratory)



Vertical Pressure Steam Sterilizer

Technical Data


MODEL	LS-35HG	LS-50HG	LS-75HG	LS-100HG
Technical Data				
Chamber volume	35L φ 318 × 450mm	50L φ 340 × 550mm	75L φ 400 × 600mm	100L φ 440 × 650mm
Working pressure	0.22MPa			
Working temperature	134°C			
Max. working pressure	0.23MPa			
Heat average	≤ ± 1°C			
Timer range	0~99min or 0~99hour59min			
Temperature adjust range	0~134°C			
Power	1 × 2.5KW/AC220V.50Hz	2 × 1.5KW/AC220V.50Hz	3 × 1.5KW/AC220V.50Hz	
Overall dimension(mm)	470 × 450 × 1020	710 × 465 × 1150	750 × 510 × 1250	750 × 510 × 1250
Transport dimension(mm)	710 × 490 × 1130	740 × 550 × 1280	780 × 600 × 1380	780 × 600 × 1380
G.W/N.W	104kg / 85kg	122kg / 99kg	146kg / 120kg	148kg / 122kg

High quality of an international norm

As a provincial high-tech enterprise, Binjiang Medical, who relies on the accurate positioning and professional spirit of being brave enough to fight, creates a miracle of development speed in the industry and wins satisfaction and trust from customers.



Professional Accreditation

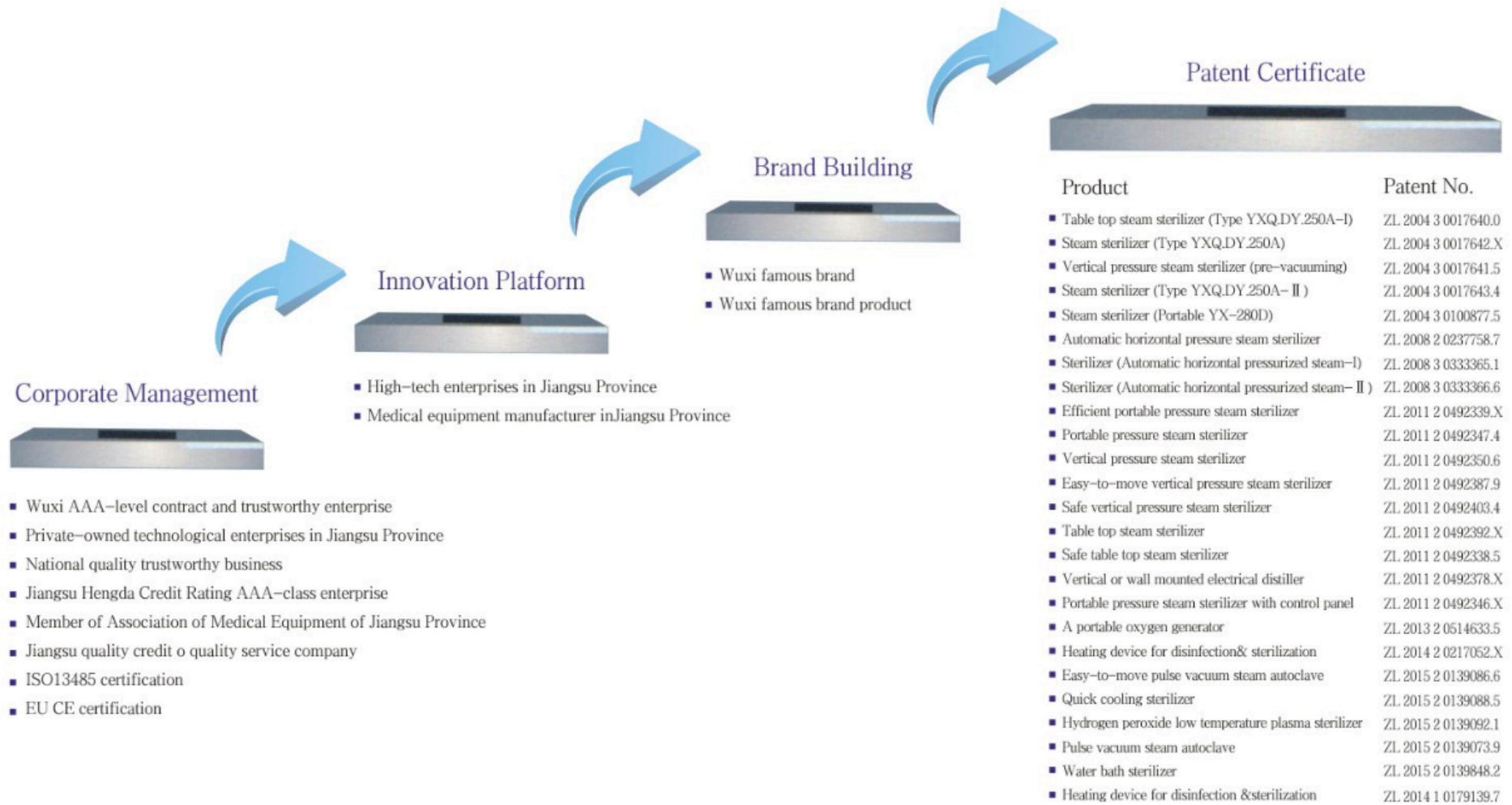

JIBI MED[®]

Jiangyin Binjiang Medical Equipment Co., Ltd.

Add: No. 38, Yungu Road, Changshou Community,
Zhouzhuang Town, Jiangyin City, Jiangsu Province

National toll free hotline: **400-700-6180**

Tel: +86-510-86270699 86296768

Fax: +86-510-86270599

P.C.: 214400

Http://www.chengyi-cn.com

E-mail: postmaster@chengyi-cn.com

jibimed@vip.163.com

