

JIBI MED[®]

DEVOTE OURSELVES TO THE PROFESSIONAL MANUFACTURING OF THE STERILIZERS

| Company Profile

Jiangyin Binjiang Medical Equipment Co., Ltd(JIBIMED) is a professional high-tech enterprise, incorporating research, manufacturing, distribution of sterilizing equipment and relevant medical products, located in Jiangyin city of Jiangsu Province, which is a remarkable place with picturesque scenery. Jiangyin is the birthplace of outstanding people, and ranked first on the economy among all the Chinese cities.

With the foundation concept of "technology reflects our wisdom, quality reflects our dignity". JIBIMED has gathered many top professional experts and famous consultants who involving themselves in the research and development of the sterile products after 10 years development. We have formed strong power to follow the tip technology at home and abroad, has become one of the largest and most professional manufacturers with continuous growth of export for the past years.

At present, JIBIMED has two producing zones, covering the area of 40 acres and the building area of 30000 m², we have more than ten full producing line for the vertical pressure steam sterilizer, portable pressure steam sterilizer, table top steam sterilizer, horizontal pressure steam sterilizer, low temperature plasma sterilizer, ultraviolet lamp trolley and intelligent drying oven, reaches 50000 sets for annual production, which are not only covered 70% of the Chinese hospitals, but also delivered to Europe, Africa, East-South Asia, Middle East over 50 countries, we have provided our goods and service to amount of tens of thousands of hospitals all over the world.

All the products from JIBIMED are strictly treated according to quality control system, such as CE, ISO9001, ISO13485. To supply the "zero-defect" units with comprehensive, considerate and sincere service is our final target goal.

The customers satisfaction and the quality expectation is our unswerving pursuit. Advanced technology, scientific management and excellent service is the guarantee of the continuously improving for the quality. Product safety, regulatory compliance and environmental protection is our social responsibility.

JIBIMED has win highly appreciation among domestic and oversea hospitals of various levels in the sterilizing field with innovated products, qualified service, and efficient management, has been awarded the "National Quality Trustworthy Enterprise", "Jiangsu Integrity and Outstanding Enterprises", ect .

For the honor and the commitment , we strive to do better !



Vertical Pulse Vacuum Steam Sterilizer

Pulse vacuum automatic



Use

The vertical pressure steam sterilizers are assembled with heating system, micro-computer controlled system, over heat and over pressure protection system, which are reliable for the sterilizing effects, convenient of operation and energy conservation. They are ideal equipments for clinics, scientific research institutes and other organizations to sterilize surgical instruments, fabrics, glasses, culture media etc.

Characteristics

- ◆ Microcomputer controlled, digital display of working statues, touch type key.
- ◆ It dynamically shows the temperature, pressure time ,vacuum rate and the cycle status.
- ◆ With fully automatically route of the water inlet, pulse vacuum, heating, sterilizing, draining, exhausting and vacuum drying.
- ◆ 3 times pulse vacuum, the vacuum rate would reach -0.081MPa after that, which ensure the completely extracting of inner air from the dressing and A,B type hollowed appliance.
- ◆ With 4 fixed programs for the appliance, medical dressing, rubber and fluid, the user also could set the program as they need randomly.
- ◆ Equipped with the PT/TT testing connector, easy for running the testing.
- ◆ Multiple pulse vacuum drying after sterilizing, which ensure the good sterilizing result.
- ◆ Equipped with 0.22μm bacteria filter, it ensure the air which enter the inner vacuumed chamber is sterile, avoid the second time contamination.
- ◆ Hand wheel type for quick-open the door. Self-inflating type seal.
- ◆ Adopted with immersed heating tube, with high heating efficiency.
- ◆ Steam-water inner circulation system: no steam discharge, and the environment for the sterilizing will be clean and dry.
- ◆ Automatically shut off with beep reminding after sterilization.
- ◆ Fully high quality stainless, equipped with double stainless steel baskets.

Technical Data

MODEL	LS-50HV	LS-75HV	LS-100HV
Technical Data			
Chamber volume	50L φ 340 × 550mm	75L φ 400 × 600mm	100L φ 440 × 650mm
Working pressure		0.22MPa	
Working temperature		134°C	
Max. working pressure		0.23MPa	
Heat average		≤ ± 1°C	
Timer range		0 ~ 99min	
Temperature adjust range		0 ~ 134°C	
Power	(3.2+1.6)KW/AC220V.50Hz	(3.2+1.8)KW/AC220V.50Hz	(3.2+1.9)KW/AC220V.50Hz
Overall dimension(mm)	660 × 460 × 1160	700 × 510 × 1250	700 × 510 × 1250
Transport dimension(mm)	740 × 550 × 1270	780 × 600 × 1380	780 × 600 × 1380
G.W/N.W	137kg / 114kg	168kg / 142kg	172kg / 146kg


Note: Equipped with separate steam generator, and printer is optional.

High quality of an international norm

As a provincial high-tech enterprise, Binjiang Medical, who relies on the accurate positioning and professional spirit of being brave enough to fight, creates a miracle of development speed in the industry and wins satisfaction and trust from customers.



Professional Accreditation

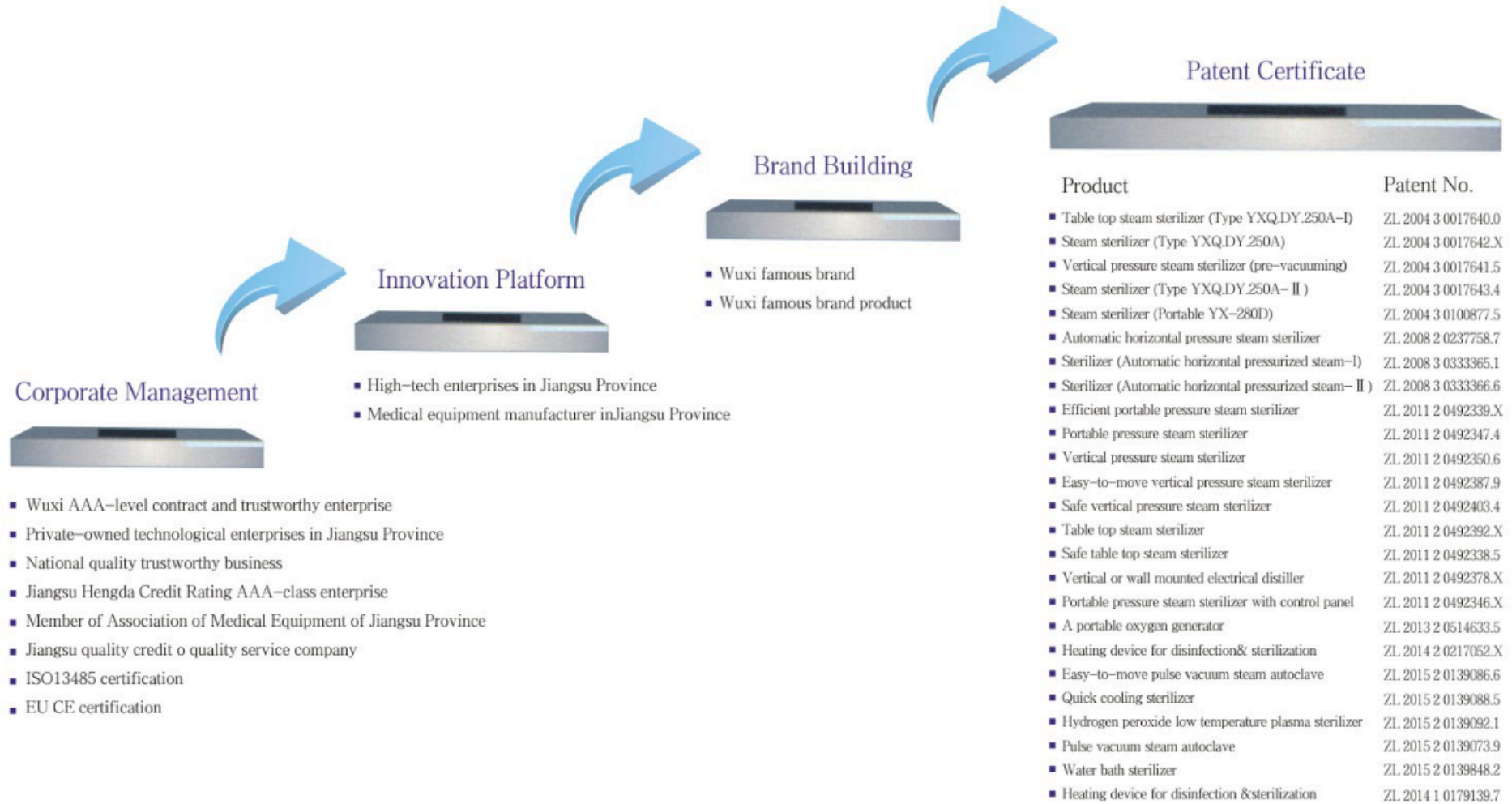

JIBI MED®

Jiangyin Binjiang Medical Equipment Co., Ltd.

Add: No. 38, Yungu Road, Changshou Community,
Zhouzhuang Town, Jiangyin City, Jiangsu Province

National toll free hotline: **400-700-6180**

Tel: +86-510-86270699 86296768

Fax: +86-510-86270599

P.C.: 214400

Http://www.chengyi-cn.com

E-mail: postmaster@chengyi-cn.com

jibimed@vip.163.com

