


Table Top Pulse Vacuum Steam Sterilizer

TM-XD12/20/24DV

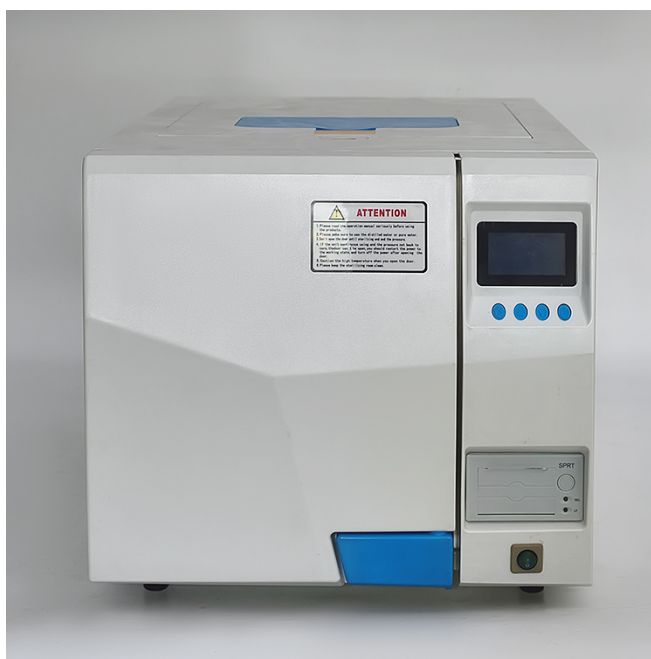

Use

The table type sterilizers are safely auto-controlled, which are designed for the dental clinics, hospitals, school lab, scientific research institutions and laboratory. And mostly suitable for sterilizing of surgical, dental and eye instruments, glassware medicine, culture medium and biological dressing, food, etc.

Characteristics

- Standard Class B with three times vacuum and drying, the remaining temperature of sterilized instruments less than 0.2%.
- LCD display working status, touch keys.
- The vacuum data can reach -0.8bar, it fits various instrument, including package, unpackaged, solid, hollow, multi apertures and inset pipes. It makes sure the instrument can be sterilized completely.
- Fully computer controlled LCD display bright show. The interface is modularized panel for easy operation.
- It is installed with BOWIE&DICK which measures the penetration of water steam.
- It has vacuum test which can make sure the ability of vacuum testing. Fast and independent steam generator.
- Built-in mini printer which can record the sterilizing information.

- Equipped with Safety valve, Safety lock door system, Pressure lock system, Pressure or temperature overload protection, Alarm system.
- It inactivates the HIV, HBV, BSE, and bacillus efficiently.
- This product is designed for open water tank, water purification, sewage tank separation, convenient and intuitive.

Technical Data

MODEL	TM-XD12DV	TM-XD20DV	TM-XD24DV
Technical Data			
Sterilizing chamber volume	12L φ200×370mm	18L φ250×360mm	24L φ250×470mm
Working pressure	0.22Mpa		
Working temperature	134℃		
Temperature adjust range	105-134℃		
Timer range	0-99min		
Heat average	≤±1℃		
Power	AC110/220V.60Hz/1.8KW	AC110/220V.60Hz/2KW	AC110/220V.60Hz/2KW
Dimension (mm)	560×470×400	560×470×400	675×470×400
Exerior package dimension (mm)	655×560×495	655×560×495	765×560×495
G.W/N.W	48/43Kg	52/47Kg	58/53Kg

Details

