


**BIOBASE<sup>®</sup>**

## **Disinfection & Sterilization Solution**

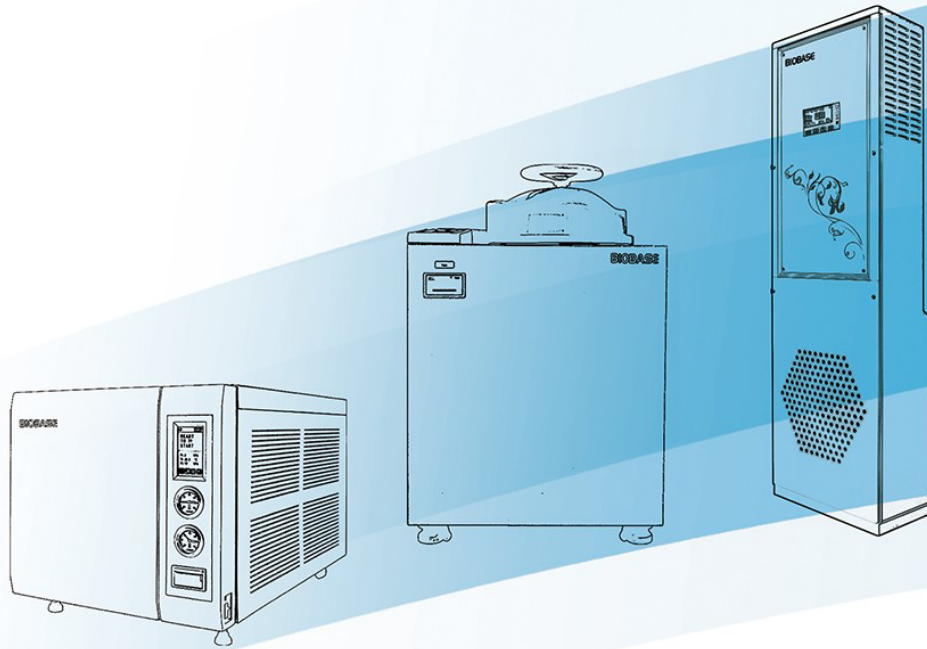

Version 2021

# About Us

Founded in 1999, having over 2000 employees (till August, 2018), BIOBASE Group is a new high-tech enterprise focusing on development, production and management of scientific equipment, In Vitro Diagnosis (IVD) instrument and reagent.

BIOBASE Group is specialized in products of 8 areas including medical diagnosis, biosafety protection, disinfection and sterilization solution, water purification system, infant care products, cold chain products, software products, clean room project.

We have full production line for table top steam sterilizer, vertical pressure steam sterilizer, portable pressure steam sterilizer, horizontal pressure steam sterilizer, hot air sterilizer, gas sterilizer, glass bead sterilizer, infrared sterilizer and automatic glassware washer. And every year we have many new products that come into the market.

BIOBASE Group has certificates including ISO9001, ISO13485, ISO14001, SGS, CE, etc for various products.

BIOBASE has set up branch offices in 24 provinces in China, and has also set up warehouses in the United States, India, UAE and Jordan. We now have long-term business partners in more than 130 countries.

BIOBASE Group always adheres to the mission of "service for health for the future" and promotes company's culture vigorously. Providing with high quality products and good service, BIOBASE would like to establish win-win business cooperation with distributors all over the world.



ISO9001



ISO13485



ISO14001



CE



## Ozone Air Sterilizer



BKCY-Y-20G



BKCY-B-20G

### Applications:

Ozone Air Sterilizer uses the strong oxidizing properties of ozone gas to make pathogenic microorganisms lose their vitality, eventually dissolve and die. Ozone disinfection has a good disinfection effect on common bacterial spores, viruses, fungi, etc., but because of its toxicity, it cannot coexist with humans and machines.

It is suitable for static disinfection of indoor air in families, schools, libraries, blood banks, nursing homes, museums, archives, restaurants, public places, etc.

### Features:

- \* Built-in high-concentration heat-dissipating generator, high sterilization efficiency.
- \* High-brightness digital tube display is used to visually display the working status of the sterilizer.
- \* The appearance is small and light, and it is made of metal material through modern moisture-proof technology, so it does not hide dirt.
- \* Infrared remote control operation mode is adopted, which is convenient for users to operate and safer.

### Technical Parameters:

Model	BKCY-Y-20G	BKCY-B-20G
Type	Mobile	Wall mounted
Disinfection Type	Ozone	
Ozone Generation Capacity	20g/h	
Suitable Disinfection Volume	60m <sup>3</sup>	
Power	≤130w	
Power Supply	Standard: 220V 50Hz, optional: 220V 60Hz, 110V 50/60Hz	
External Size	198*206*485mm	
Net Weight	7.2kg	6.8kg
Packing Size	274*286*567mm	
Gross Weight	8.8kg	8.3kg

## **BIOBASE GROUP**

Add: No.9 Gangxing Road, High-tech Zone,  
Jinan City, Shandong Province, China  
Tel: +86-531-81219801/03 Fax: +86-531-81219804  
E-mail: [export@biobase.com](mailto:export@biobase.com)  
Website: [www.biobase.cc](http://www.biobase.cc) / [www.biobase.com](http://www.biobase.com)



BIOBASE CHINA



BIOBASE CHINA



BIOBASE CHINA