


BIOBASE[®]

Disinfection & Sterilization Solution

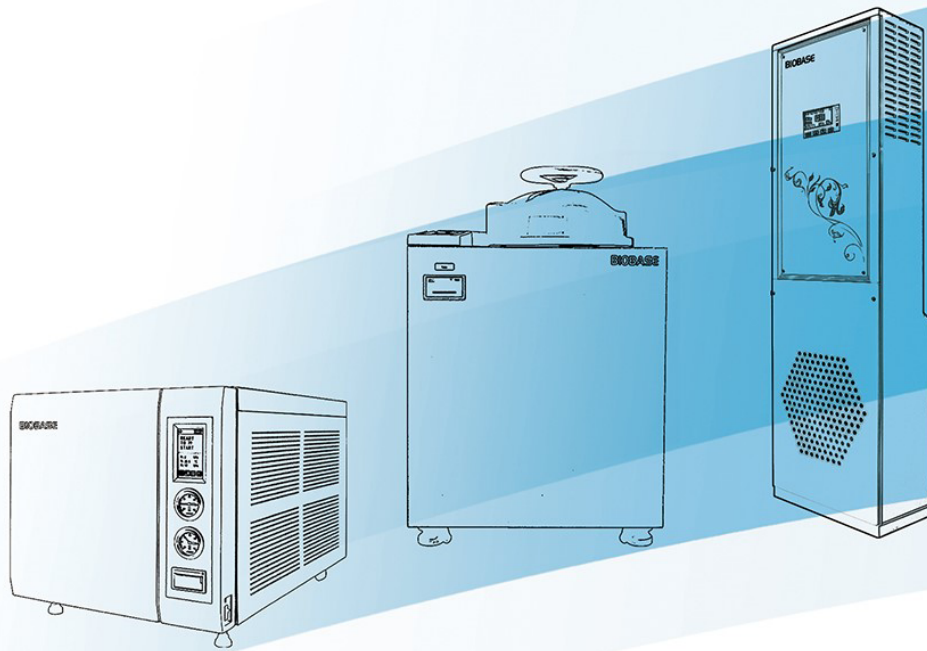

Version 2021

About Us

Founded in 1999, having over 2000 employees (till August, 2018), BIOBASE Group is a new high-tech enterprise focusing on development, production and management of scientific equipment, In Vitro Diagnosis (IVD) instrument and reagent.

BIOBASE Group is specialized in products of 8 areas including medical diagnosis, biosafety protection, disinfection and sterilization solution, water purification system, infant care products, cold chain products, software products, clean room project.

We have full production line for table top steam sterilizer, vertical pressure steam sterilizer, portable pressure steam sterilizer, horizontal pressure steam sterilizer, hot air sterilizer, gas sterilizer, glass bead sterilizer, infrared sterilizer and automatic glassware washer. And every year we have many new products that come into the market.

BIOBASE Group has certificates including ISO9001, ISO13485, ISO14001, SGS, CE, etc for various products.

BIOBASE has set up branch offices in 24 provinces in China, and has also set up warehouses in the United States, India, UAE and Jordan. We now have long-term business partners in more than 130 countries.

BIOBASE Group always adheres to the mission of "service for health for the future" and promotes company's culture vigorously. Providing with high quality products and good service, BIOBASE would like to establish win-win business cooperation with distributors all over the world.



ISO9001



ISO13485



ISO14001



CE



H₂O₂ Nebulizer Sterilizer



Applications:

It is suitable for biosafety laboratories, infection wards, schools, hospitals, food enterprises, scientific research institutes and other places that require regular disinfection of air and hard surfaces.



Features:

- * It can control the spraying of the disinfectant in the form of nano-micron droplets, reduce the condensed liquid, and ensure the penetration and disinfection effect on complicated and difficult-to-reach places.
- * Full-color touch screen, easy and clear operation, high user experience comfort.
- * It can achieve a killing effect of 6 logs of Bacillus subtilis black mutated spores.
- * It is less corrosive, non-toxic, and has a wide range of applications, which can replace traditional space disinfection methods.
- * Flexible to use and easy to operate, small in size and easy to move, suitable for various units that have requirements for microbial limits. It has an automatic calculation function. Just input the size of the room to be disinfected, and the required amount of disinfectant can be automatically calculated.
- * With the function of timing appointment, this product can make an appointment in advance. After the timing appointment is made, the equipment can automatically start disinfection.

Technical Parameters:

Model	BKX-G100
Liquid Capacity	4L
Suitable Space	100m ³
Power	1100W
Power Supply	220V 50Hz
External Size	530*530*760mm
Net Weight	48kg
Packing Size	640*640*885mm
Gross Weight	55kg

BIOBASE GROUP

Add: No.9 Gangxing Road, High-tech Zone,
Jinan City, Shandong Province, China
Tel: +86-531-81219801/03 Fax: +86-531-81219804
E-mail: export@biobase.com
Website: www.biobase.cc / www.biobase.com



BIOBASE CHINA



BIOBASE CHINA



BIOBASE CHINA