

JIBI MED[®]

DEVOTE OURSELVES TO THE PROFESSIONAL MANUFACTURING OF THE STERILIZERS

| Company Profile

Jiangyin Binjiang Medical Equipment Co., Ltd(JIBIMED) is a professional high-tech enterprise, incorporating research, manufacturing, distribution of sterilizing equipment and relevant medical products, located in Jiangyin city of Jiangsu Province, which is a remarkable place with picturesque scenery. Jiangyin is the birthplace of outstanding people, and ranked first on the economy among all the Chinese cities.

With the foundation concept of "technology reflects our wisdom, quality reflects our dignity". JIBIMED has gathered many top professional experts and famous consultants who involving themselves in the research and development of the sterile products after 10 years development. We have formed strong power to follow the tip technology at home and abroad, has become one of the largest and most professional manufacturers with continuous growth of export for the past years.

At present, JIBIMED has two producing zones, covering the area of 40 acres and the building area of 30000 m², we have more than ten full producing line for the vertical pressure steam sterilizer, portable pressure steam sterilizer, table top steam sterilizer, horizontal pressure steam sterilizer, low temperature plasma sterilizer, ultraviolet lamp trolley and intelligent drying oven, reaches 50000 sets for annual production, which are not only covered 70% of the Chinese hospitals, but also delivered to Europe, Africa, East-South Asia, Middle East over 50 countries, we have provided our goods and service to amount of tens of thousands of hospitals all over the world.

All the products from JIBIMED are strictly treated according to quality control system, such as CE, ISO9001, ISO13485. To supply the "zero-defect" units with comprehensive, considerate and sincere service is our final target goal.

The customers satisfaction and the quality expectation is our unswerving pursuit. Advanced technology, scientific management and excellent service is the guarantee of the continuously improving for the quality. Product safety, regulatory compliance and environmental protection is our social responsibility.

JIBIMED has win highly appreciation among domestic and oversea hospitals of various levels in the sterilizing field with innovated products, qualified service, and efficient management, has been awarded the "National Quality Trustworthy Enterprise", "Jiangsu Integrity and Outstanding Enterprises", ect .

For the honor and the commitment , we strive to do better !



Plasma Air Sterilizer (Mobile Type)



Use

JIBIMED series plasma air sterilizer mainly uses plasma as sterilization factor, which can be used to sterilize indoor air, kill bacteria in air, and can be used in human places. It is suitable for air disinfection in II, III, IV environment of general operating room, delivery room, baby room, premature delivery room, supply room, burn ward and so on, as well as indoor air disinfection in pharmaceutical, food, school, public places and other industries.

Characteristics

- ◆ Filtration and dust removal: the physical filtration method in air laminar cleaning technology is used to remove dust and bacteria from the air.
- ◆ Safe and reliable: this machine is dynamic air disinfection equipment, man-machine coexistence, harmless to people, no corrosion to equipment, but also can be used for static air disinfection in unmanned places.
- ◆ Convenient to use: All the machines have the functions of microcomputer program control, far infrared remote control, automatic monitoring of plasma fault, wind speed regulation, automatic accumulation of working time and so on.

Technical Data


	MODEL	
Technical Data	BJ-100Y	BJ-120Y
Index applicable volume(m ³) ≤	100	120
Suggested time(min) with high wind speed	90	90
Circulating air volume(m ³ h) ≥	800	960
Input power(VA) ≤	150	155
Fuse(A)	2	
Noise db(A) ≤	45	
Plasma density	4.08 × 10 ¹⁶ /m ³ ~1.37 × 10 ¹⁷ /m ³	
Ozone content in working environment(mg/m ³)	0.02	
Disinfection efficiency	Within 120 minutes, the extinction rate of the simulated experiment of Staphylococcus billycans is ≥99.9%, and the extinction rate of the natural bacteria in the air in the field experiment is ≥90%.	
Installation size(cm)	36 × 35 × 97	42 × 40 × 101
Installation mode	Mobile type	

High quality of an international norm

As a provincial high-tech enterprise, Binjiang Medical, who relies on the accurate positioning and professional spirit of being brave enough to fight, creates a miracle of development speed in the industry and wins satisfaction and trust from customers.



Professional Accreditation

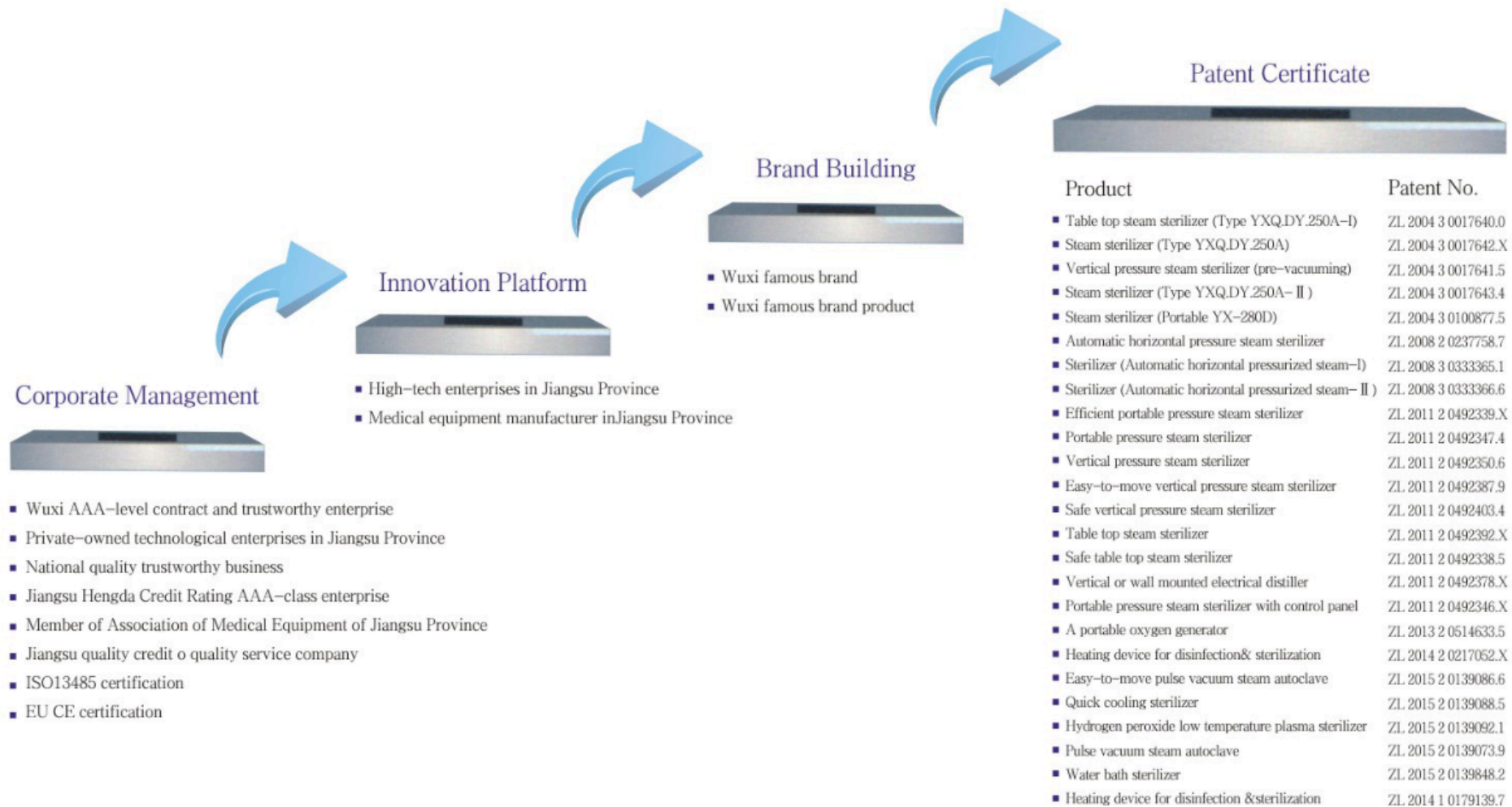

JIBI MED®

Jiangyin Binjiang Medical Equipment Co., Ltd.

Add: No. 38, Yungu Road, Changshou Community,
Zhouzhuang Town, Jiangyin City, Jiangsu Province

National toll free hotline: **400-700-6180**

Tel: +86-510-86270699 86296768

Fax: +86-510-86270599

P.C.: 214400

Http://www.chengyi-cn.com

E-mail: postmaster@chengyi-cn.com

jibimed@vip.163.com

