



# SKB041-3S

## HYDRAULIC PATIENT TROLLEY

Registration code: SKM-G




**SAIKANG** Medical

Healthier  
your life

# SKB041-3S

## FEATURES AND HIGHLIGHTS

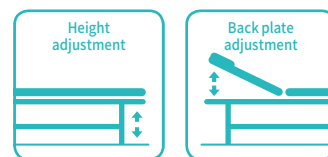

### Technical parameters

▪ External size (LxW):	2150×790mm
▪ Table height adjustment:	635-895mm
▪ Back plate adjustment:	0-60° (±10°)

### Technical configuration

☑ Hydraulic Pump	1 pc
☑ Air Spring	1 pc
☑ Stainless steel guardrail	1 set
☑ Push handle	1 set
☑ Wall bumper wheel	4 pcs
☑ Foot pedal	2 pcs
☑ 5" Castors with brake	4 pcs
☑ Oxygen cylinder rack	1 pc
☑ Debris basin	1 pc
☑ IV pole	1 pc
☑ Mattress with thickness of 60mm	1 pc



**Metal material**

Baosteel from Fortune Global 500  
Panasonic Robotic ensure 360° full  
smooth welding.

**Painting**

11 process epoxy painting, ASTM  
testing anti-bacterial, paint thickness  
0.12 mm, brightness 60°. Painting can  
resist 50 kg impact.

**Multifunctional patient cart**

Multi-functional trolley, which can  
be separated into two parts: stretcher  
and trolley.  
The stretcher with four wheels can be  
used and moved individually.

**Bumper wheel**

Bumper wheel detachable, max.  
protect the bed and patient during  
transportation.

**Guardrail**

Foldable stainless steel guardrail, no  
rust, easy to clean.

**Hydraulic pumps**

American brand hydraulic pump,  
providing strong support for the  
functions of height adjustment.

**Casters**

Four covered casters of 5 inch  
diameter, all with brakes.

**IV pole**

Foldable telescopic stainless steel IV  
pole, bearing 15kg.

**Oxygen bottle holder**

Oxygen bottle holder, load capacity:  
15KG.

**Stainless steel basin**

Equipped with stainless steel sewage  
basin and other various accessories.

**Mattress**

Waterproof mattress, convenient for  
transfer patient.





CONSTANT  
INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE  
BRANCHES

#### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

#### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

#### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

#### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

#### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SAIKANG** Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎ +86 18021231901 · 🕒 24/7  
✉ sales@chinasaikang.com  
🌐 www.saikangmedical.com



Company video



General catalogue